

**OLYMPUS**<sup>®</sup>

Your Vision, Our Future



EVIS EXERA GASTROINTESTINAL VIDEOSCOPE  
**OLYMPUS GIF TYPE XP160**

# SUPER-SLIM INSERTION TUBE

With its ultra-narrow 5.9 mm diameter insertion tube, this SlimSIGHT™ scope is easier to insert into the upper digestive tract while its convenient distal-end irrigation nozzle ensures clear images

**SlimSIGHT**<sup>™</sup>




## Proven Features for Powerful Performance

- Extra-narrow insertion tube measuring a mere 5.9 mm across.
- Ultra-slim design ideal for patients with strictures and narrow pharynges or esophagi, as well as for routine examinations.
- Continuously clear view ensured by distal-end irrigation nozzle that washes away any material that adheres to the lens during examination.
- 2.0 mm diameter channel incorporated for suctioning and biopsies.
- 120° field of view and 4-way angulation (180° up, 90° down and 100° right/left) enable complete examination of the upper digestive tract.
- Ergonomically designed grip enhances scope maneuverability while redesigned easy-to-access knobs and user-selectable switches improve operability.
- Superior image quality is delivered even in the expanded full height mode when combined with the dedicated CV-160 video system center.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.



### Specifications

Optical System	Field of view	120°
	Direction of view	0° Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	5.9 mm
	Insertion Tube	Outer diameter
Bending Section	Angulation range	Up 180°, Down 90°, Right 100°, Left 100°
Working Length		1030 mm
Total Length		1345 mm
Instrument Channel	Inner diameter	2 mm
	Minimum visible distance	3 mm from the distal end
	Endo-Therapy accessory entrance/exit position in field of view	

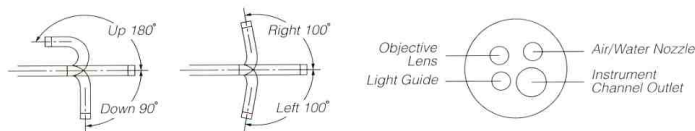

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

**OLYMPUS CORPORATION**  
2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan  
**OLYMPUS EUROPA GMBH**  
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany  
**OLYMPUS AMERICA INC.**  
2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.  
**OLYMPUS LATIN AMERICA, INC.**  
6100 Blue Lagoon Drive, Suite 390 Miami, Florida 33126-2087, U.S.A  
**KEYMED LTD.**  
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

**OLYMPUS SINGAPORE PTE LTD.**  
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373  
**OLYMPUS HONG KONG AND CHINA LIMITED**  
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
**OLYMPUS BEIJING REPRESENTATIVE OFFICE**  
Room 1401, East Tower, Gong Yuan No6 Royal Palace, No6 Gon YuanXijie, JianGuoMenNei, DongCheng District, Beijing, 100005, China  
**OLYMPUS MOSCOW LIMITED LIABILITY COMPANY**  
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia  
**OLYMPUS AUSTRALIA PTY LTD.**  
31 Gilby Road, Mount Waverley, VIC., 3149, Australia